

FOCUSLIFT



1

Lighten
wrinkles

2

lifting &
tightening

3

promote new
students

Operator Interface




Handle



Handle



Key switch

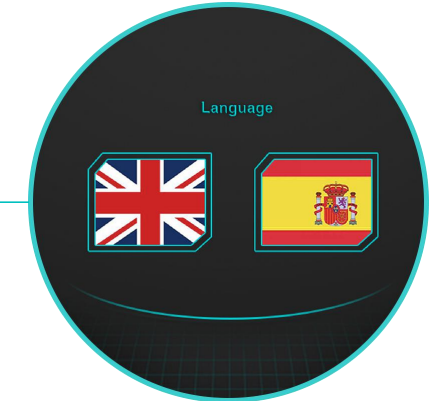

Touch screen

Operator Interface



Content

Latest launched High intensity focused Radio Frequency skin lifting system

New powerful focused RF, Interval heating depth focus.


Radar positioning and navigation technology. Stratify meticulously, layered heating, 3 layers tighten.

Intelligent recognition: Intelligent identification the collagen and fat, no damage to other normal tissues.

RF with vacuum and cooling, more comfortable and better results




Treatment head

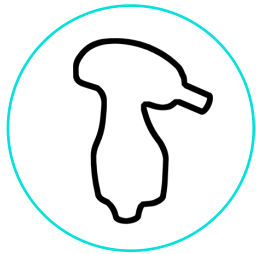

Mono polar RF

Mono polar RF treatment head with contact cooling for face treatment

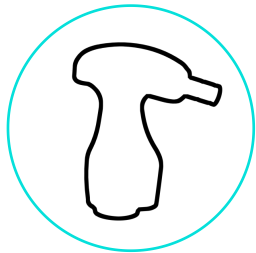


Mutiploar RF

Mutiploar RF treatment head with pulsed vacuum and contact cooling for body treatment




Principle


Superficial Heat

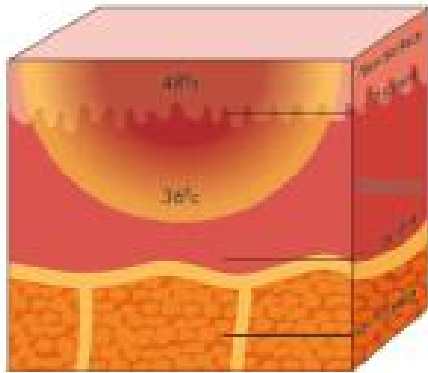

Medium Heat

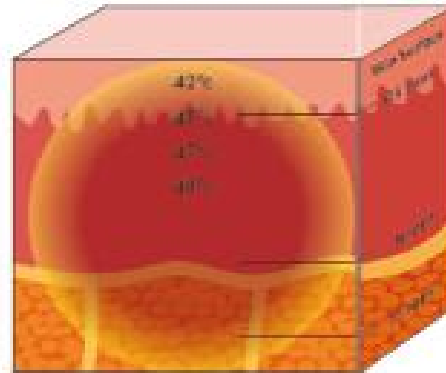


Deep Heat

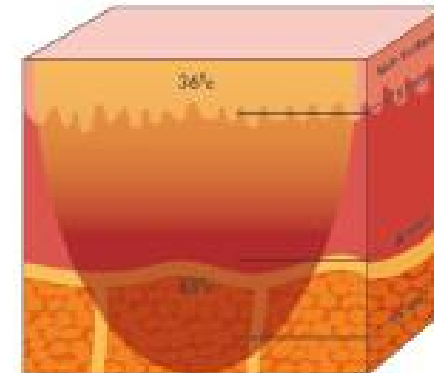

Focuslift adopts the new generation of RF firming technology, and the accurate action on epidermis layer, derma layer and SMAS layer to achieve the purpose of wrinkles removal, skin tightening, beautifying and body shaping.

Principle

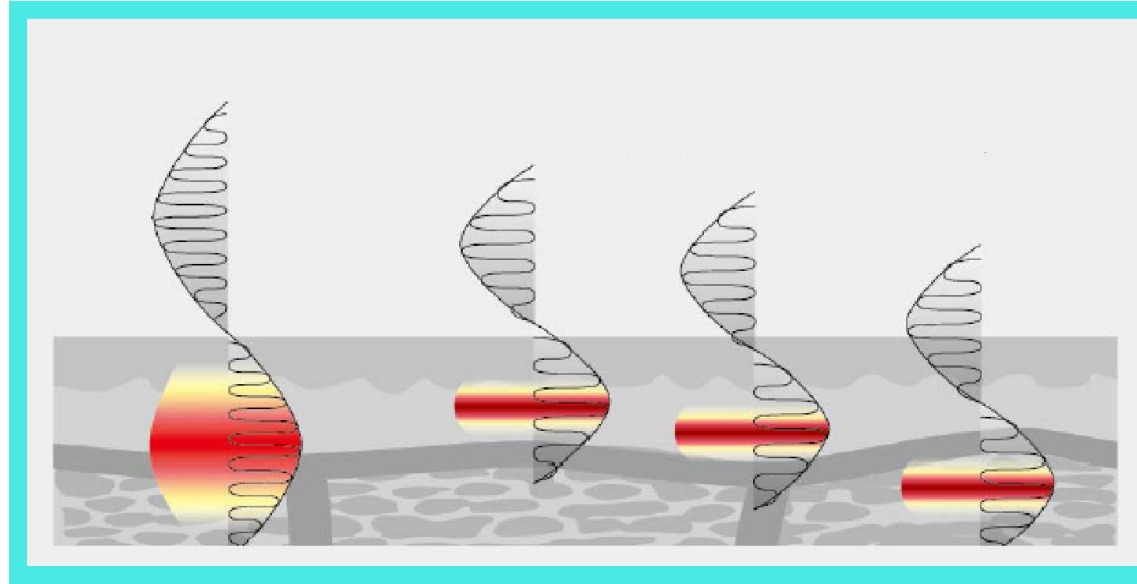

Epidermis

Dermis layer


SMAS layer

Focuslift can make the RF energy reach to the dermis accurately .Then it can promote contraction of shallow collagen and lift the skin structure immediately .


Focuslift can let the RF energy reach to the deep dermis selectively. It excites the wound healing response and release the collagen of dermis fibroblasts. Then recovering skin elastic and keeping skin tightening .

Focuslift can make the RF energy reach to the SMAS accurately. Through the thermal effect renew the elastic fiber structure and tighten the SMAS.


Contrast

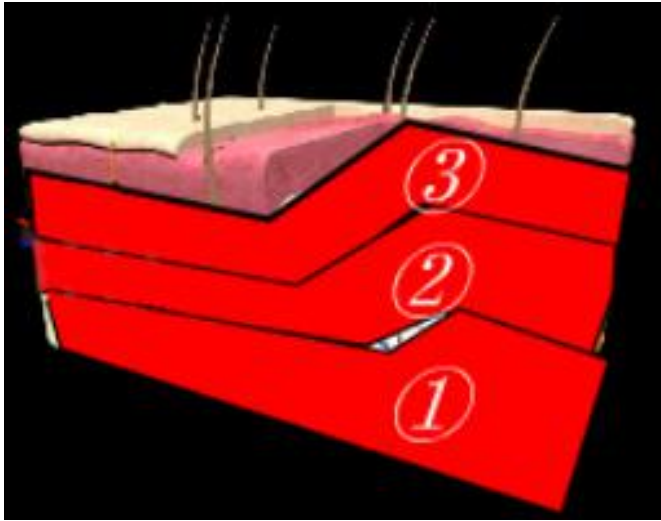

Focuslift

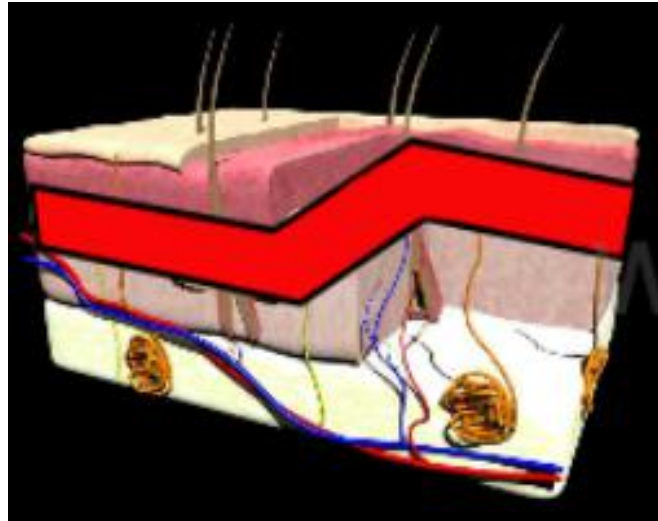

Bipolar RF

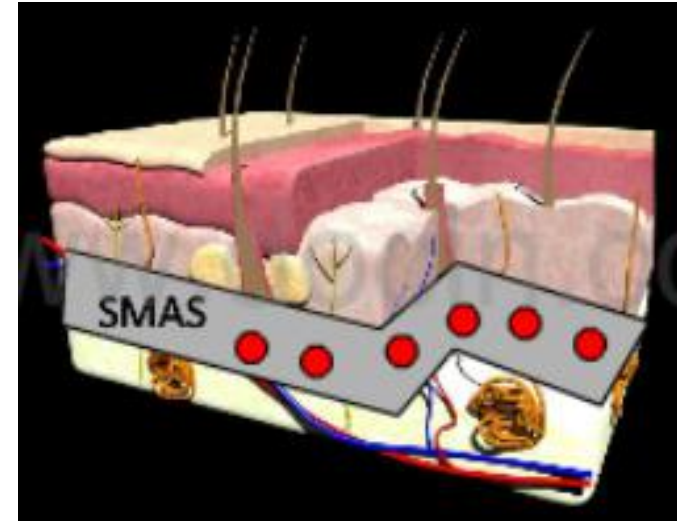

HIFU

Applications

03

Wrinkle removal : Forehead wrinkles, eye round fine lines, neck/hands/body all kinds of wrinkle.

02

Face lifting : Facial relaxation, double chin, chin brow vertical curve and brow, small eyelid tightening.

01

Body shaping : Arm tightening and shaping, back tightening and shaping, breasts shaping, waist and abdomen tightening.



Firming skin




Improve mandibular skin folds



Leaves skin smooth and shiny



Contour shaping

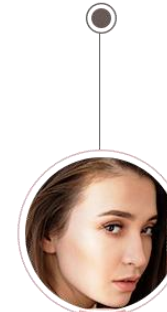

Improve skin elasticity



Eliminate forehead and around the face Wrinkles



Parameters

Product name	Foucslift skin lifting system
Power source	High intensity Focused Radio frequency
Frequency	40.68Mhz
Energy	1-10 degree
Power duration	100-1000ms, 50ms step controlling
Cooling system	Built-in water and semiconducting cooling system
Vacuum Settings	Pulsed/continued vacuum
Voltage	AC220V ±22V ,180VAAC110V+22V,180VA
Net weight	14kg
Packaging Size	460*620*600mm
Warranty	2 year for whole unit
Consumable part	No



Treatment



1



2



3



4

3 treatment with 1 month interval continually is necessary.
2-3 years continually effect will be available after 1 section
Extra treatment is necessary for long lasting better effect.

Do not wash treatment area with hot water in 3 days
Moisturizing and sun block cream is necessary after treatment