

# MACHINE DETAILS



### DETAILED DESCRIPTION



#### **4** slimming technologies

Cryolipolysis, Vacuum, RF, and Cavitation 4 slimming technologies in one makes treatment more safe and cost-effective.



#### **5** handles are available

(2 Cryo heads, 2 RF heads: one for face and one for body, 1 Cavitation head). It can do whole body treatment.



### Popular and Hot Design

Smart and high standard design of the cryo machine will help to improve the image of your salons or clinics.



#### Best Freeze -15 Celsius degree

Strong cooling system ensure the temperature reach -15degree, much lower than -10 degree in the current market.

### Super strong Cooling System

Powerful Booster water pump with unique 8L big water tank, and Japanese imported pure cooper radiator with 4 big DC fans all ensure the machine keep working continuously for 12 hours.




### **Double Air Pumps**

Equipped with 2 powerful air pumps to ensure the high performance of two handles working together.

### HANDLES

Six handles of three different sizes, either two handles could at the same time.

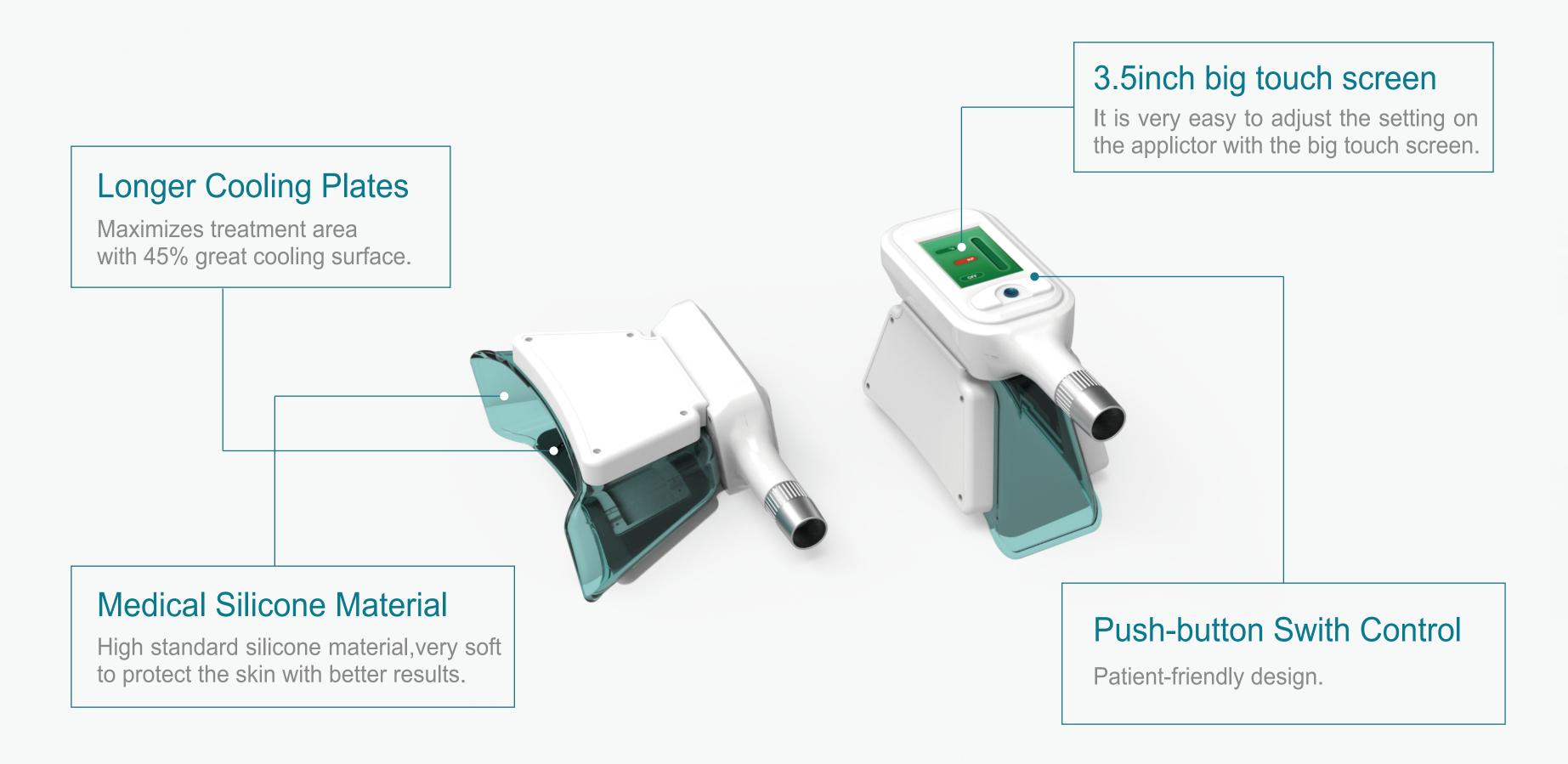


## HANDLES DETAILS





3 minutes heating before cooling to avoid frosthite

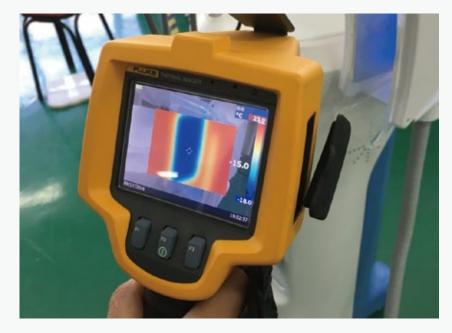

Lowest temperature -15°C



Medical use soft silicone

## TECHNICAL PARAMETERS

Model	LB-C5
Screen	12 inch touch color screen
Heads	5 pieces available (2 cryo, 2 RF, 1 cavitation)
Input Power	2000W
Cryo	Vacuum: 0kpa-100kpa
	Cooling temperature: -15-5°C
	Two cryo heads work together
RF	RF frequency: 2MHZ
	RF energy: 1- 50J/cm2
Ultrasound	Ultrasound cavitations frequency: 40K
Voltage	Ultrasound cavitations power: 1-70J/cm2
Cooling system	AC220V±10%,50Hz/ AC110V±10%,60Hz
	Wind cooling + water cooling+ semiconductor cooling


#### HOW NON-INVASIVE CRYOLIPOLYSIS WORKS

1




Some stubborn fat bulges are immune to diet and exercise

4

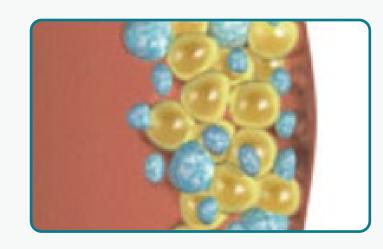


Following treatment. fat cells enter an apoptotic, death sequence and are gradually removed in the next few weeks and months by the immune system

7

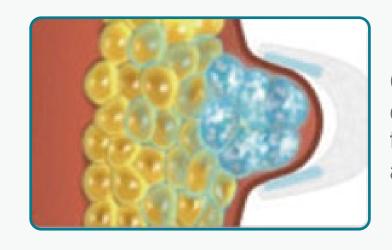


CoolSculpting targets and cools fat cells to temperatures that trigger fat cell apoptosis

5

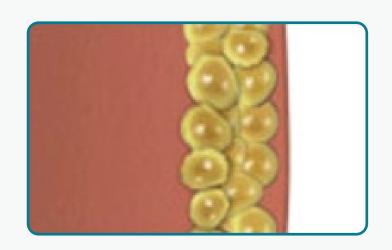

Fatlayer thickness significantly reduced

3

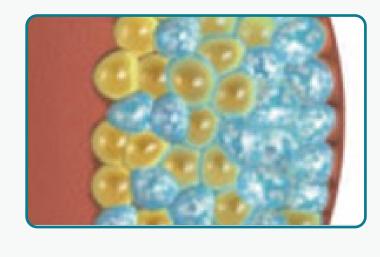


No damage to nerves or other tissues because lipids in fat crystallize at a warmer temperature than water other cell types


6




Fatlayer reduction in targeted area leads an improvement in flank appearance


## Contrast before and after

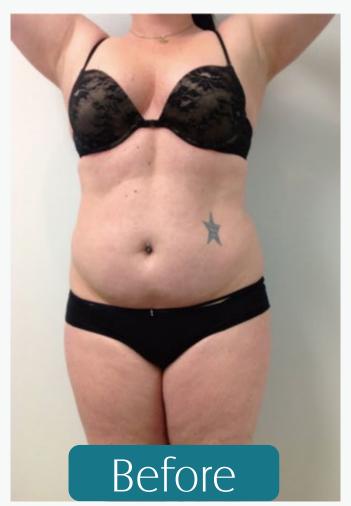

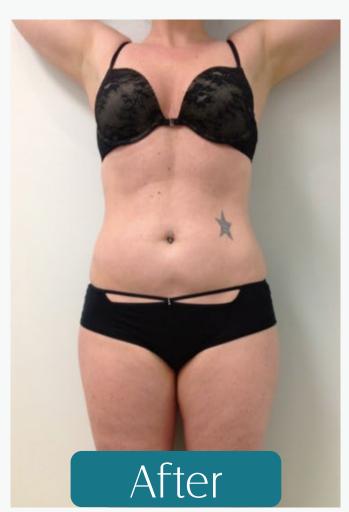


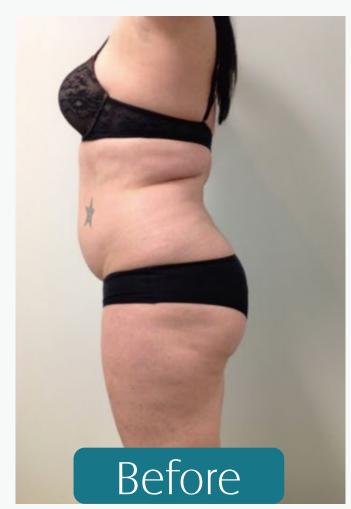



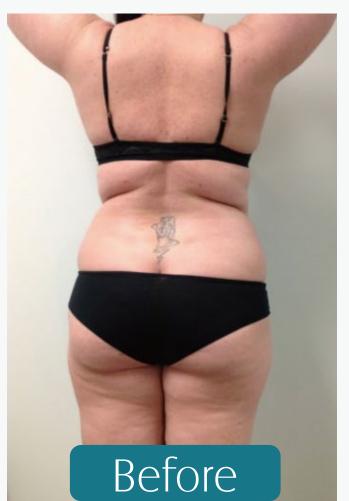



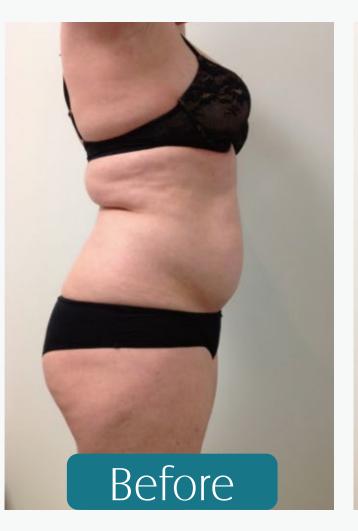




# Customer use





### Client Feedback

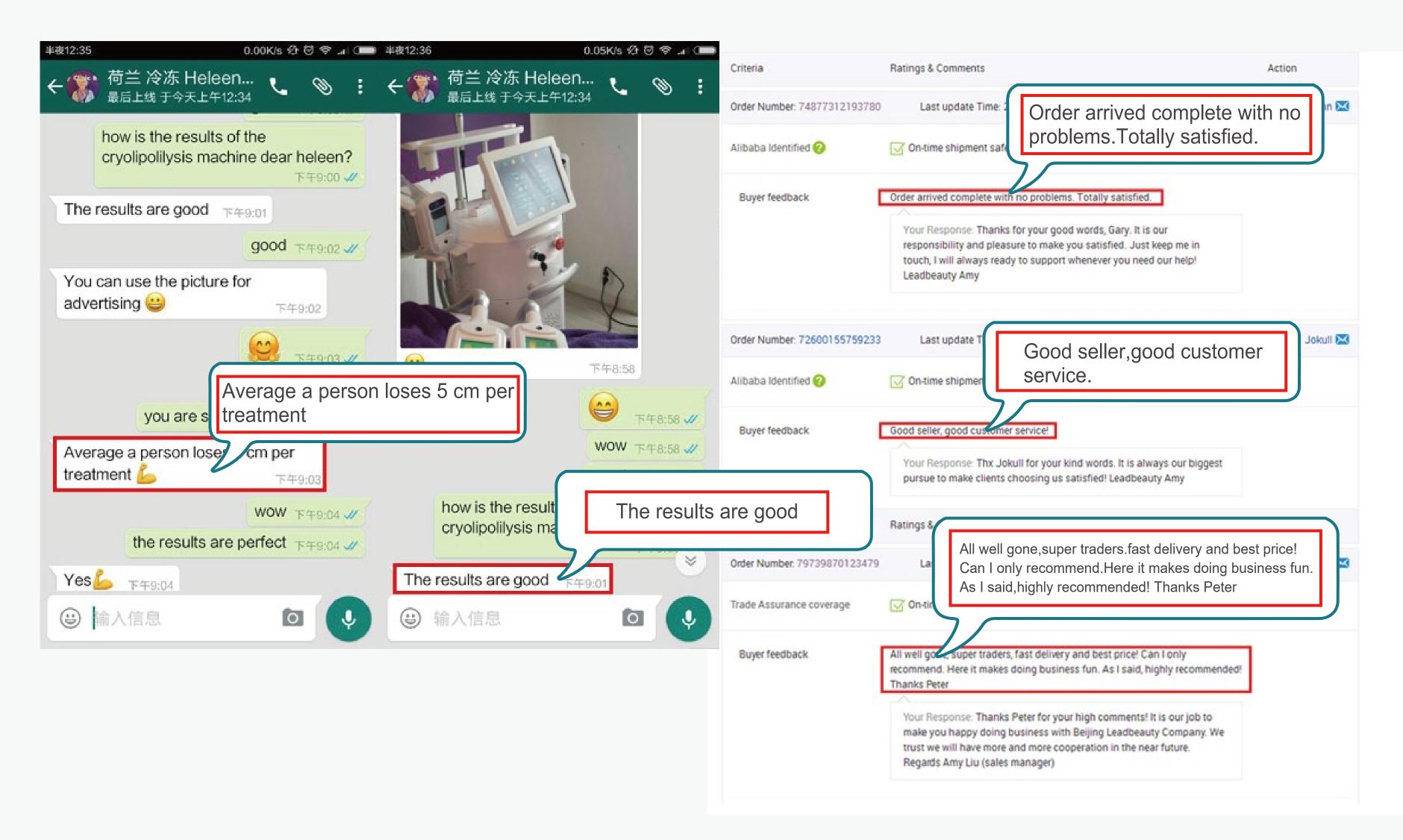

# About us Over 13 years manufacturer- Strong production capacity Beijing Leadbeauty International S&T CO. Ltd can keep you always ahead in the aesthetics & medical devices field. LEADBEAUTY has its own research & development, sales and after-sales departments; can offer the professional technology supports at the first time. We are one of the most leading & reliable manufacturers in China, have a professional team integrating with optics, machinery and electricity that keep us ahead on in this field. Beijing LEADBEAUTY International S&T CO. Ltd Tel: 86-10-60418895 Mob & whatsapp: +86 13269056123 E-mail: leadbeauty@leadbeauty.net Website: www.leadbeauty.net

Addr: West zone No.5 Mauhwa Industry Park, No.1 Caida 3rd Street, Shunyi District, Beijing, China